




BIO-SHIELD FUMONISIN 5

Cat.Number B5548/B5596

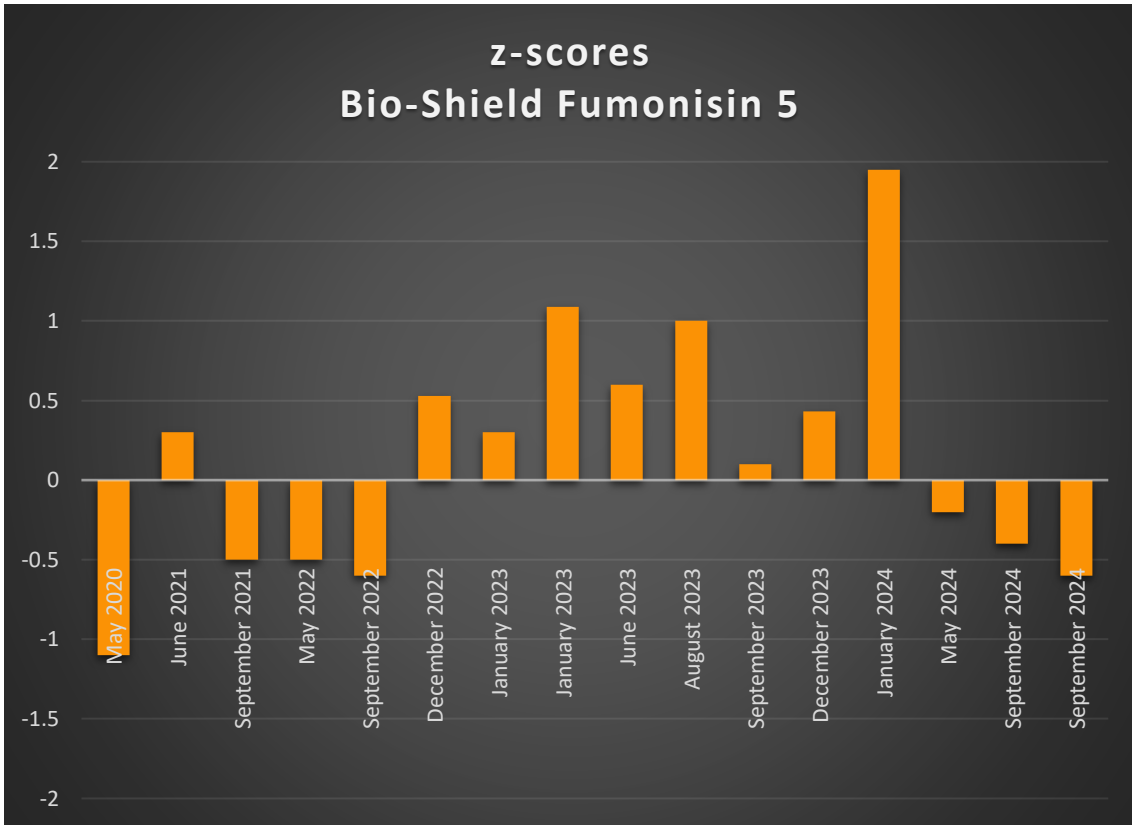

| Test | Assigned value (µg/kg) | Result (µg/kg) | Z-score |
|--|------------------------|----------------|---------|
| FAPAS Food Chemistry Proficiency Test - 04519 - Multi - Mycotoxins in Pet Dog Food (Dry, Cereal Based) - July - September 2024 | 327 | 291 | -0.6 |
| FAPAS Food Chemistry Proficiency Test - 04520 - Multi - Mycotoxins in Maize Flour - July - September 2024 | 575 | 537 | -0.4 |
| FAPAS Food Chemistry Proficiency Test - 04510 - Multi - Mycotoxins in Maize Flour - April - May 2024 | 788 | 759.07 | -0.2 |



BIO-SHIELD

| | | | |
|--|-----|--------|------|
| BIPEA 31E - Mycotoxins – Feed – Animal Feed - 14 - 5 – January 2024 | 150 | 226 | 1.95 |
| BIPEA 31B - Mycotoxins – Cereals – Ground Corn - 6 - 187 – December 2023 | 434 | 481 | 0.43 |
| FAPAS Food Chemistry Proficiency Test - 04482 – Multi - Mycotoxins in Maize Flour - July - September 2023 | 640 | 646.3 | 0.1 |
| FAPAS Food Chemistry Proficiency Test - 04481 - Mycotoxins in Pet Dog Food - July – August 2023 | 452 | 533 | 1.0 |
| FAPAS Food Chemistry Proficiency Test - 04472 - Multi - Mycotoxins in Maize Flour - April - June 2023 | 798 | 882 | 0.6 |
| BIPEA 31E - 3 - Mycotoxin Contamination in Animal feed - 13 - 3031 - January 2023 | 134 | 172 | 1.09 |
| FAPAS Food Chemistry Proficiency Test - 04461 - Mycotoxin Contamination in Animal Feed (Cereal Based) - December 2022 - January 2023 | 347 | 368.87 | 0.3 |
| BIPEA 31B - 182 - Mycotoxin Contamination in Ground corn - 18 - 0731 - December 2022 | 167 | 190 | 0.53 |
| FAPAS Food Chemistry Proficiency Test - 04448 – Multi - Mycotoxins in Maize Flour - July - September 2022 | 569 | 512.5 | -0.6 |
| FAPAS Food Chemistry Proficiency Test - 04439 – Multi - Mycotoxins in Maize Flour - April - May 2022 | 983 | 911.5 | -0.5 |



BIO-SHIELD

| | | | |
|---|-----|---------|------|
| FAPAS Food Chemistry Proficiency Test - 04418 - Multi - Mycotoxins in Maize Flour - July - September 2021 | 720 | 658.5 | -0.5 |
| FAPAS Food Chemistry Proficiency Test - 04411 - Multi - Mycotoxins in Maize Flour - April - June 2021 | 962 | 1015.13 | 0.3 |
| FAPAS Food Chemistry Proficiency Test - 04384 - Multi - Mycotoxins in Maize Flour - April - May 2020 | 752 | 612.19 | -1.1 |